

PART OF THE **DURAN**® BOTTLE SYSTEM


# **DURAN**® GL 45 BOTTLE CARRYING SYSTEM



**DURAN  
WHEATON  
KIMBLE**

Excellence in your hands

## DURAN® GL 45 BOTTLE CARRYING SYSTEM

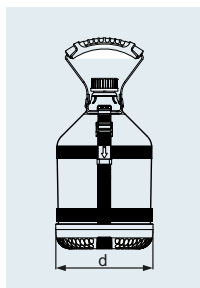

The new modular Bottle Carrying System fits standard laboratory bottles that have a GL 45 neck finish and is compatible with four bottle sizes 2, 5, 10 and 20 litres. It makes pouring liquids from the larger bottles more comfortable and safer. It is ergonomically designed, and offers health & safety benefits to the user and allows easy and safe transport of bottles around the laboratory or production facility. The collar is only removable by cutting – so attachment is 100% secure in use.

Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	Colour	Pack Unit
29 283 63 07	2 000	GL 45	136	blue	1
29 282 63 06	2 000	GL 45	136	purple	1
29 281 63 05	2 000	GL 45	136	red	1
29 280 63 04	2 000	GL 45	136	yellow	1
29 280 73 09	5 000	GL 45	194	yellow	1
29 280 86 05	10 000	GL 45	242	yellow	1
29 280 91 04	20 000	GL 45	304	yellow	1


Cat. No.	Description	Pack Unit
<b>Accessories 10 L</b>		
29 901 86 07	DURAN® Metal Dolly stainless steel for 10 L DURAN® GL 45 bottle	1
29 903 86 09	DURAN® Handle for the 10 L Metal Dolly stainless steel	1
29 280 86 54	DURAN® Silicone Bottle Base Protector for the 10 L Metal Dolly, Yellow	1
<b>Accessories 20 L</b>		
29 901 91 06	DURAN® Metal Dolly for 20 L DURAN® GL 45 bottle stainless steel	1
29 903 91 08	DURAN® Handle for the 20 L Metal Dolly stainless steel	1
29 280 91 53	DURAN® Silicone Bottle Base Protector for the 20 L Metal Dolly, Yellow	1
<b>Replacement Parts</b>		
29 284 73 04	DURAN® Replacement handle and collar for GL 45 bottle carrier 5 L	1
29 284 86 09	DURAN® Replacement handle and collar for GL 45 bottle carrier 10 or 20 L	1



## DURAN® ORIGINAL LABORATORY BOTTLE



with DIN thread

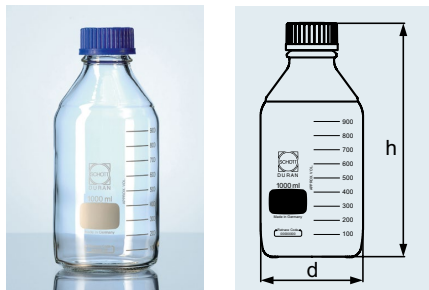

With easy-to-read scale and large labelling field for easy marking, in fired-on, highly durable white ceramic. With proven DURAN® properties. Complete with blue screw cap (PP, integral lip seal) and pouring ring (PP) for drip-free pouring and clean, safe working. Service temperature level of screw cap and pouring ring: +140 °C.


Typical applications: storage, sample preparation, transport. Autoclaving media.

Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	h (mm)	Pack Unit
<b>with screw cap and pouring ring from PP (blue)</b>					
21 801 63 57	2 000	45	136	265	10
21 801 73 53	5 000	45	182	335	1
21 801 86 58	10 000	45	227	415	1
21 801 91 57	20 000	45	288	510	1
<b>without screw cap and pouring ring</b>					
21 801 63 08	2 000	45	136	260	10
21 801 73 04	5 000	45	182	330	1
21 801 86 09	10 000	45	227	410	1
21 801 91 08	20 000	45	288	505	1

## DURAN® ORIGINAL GL 45 LABORATORY BOTTLES COMPLETE WITH HIGH TEMPERATURE CLOSURES



with DIN Thread, clear, graduated

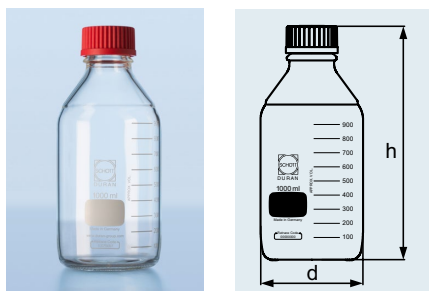

DURAN® original GL 45 Laboratory Bottles are available complete with high temperature screw caps and pouring rings. The PBT cap and ETFE pouring rings offer greater thermal and chemical resistance than the equivalent polypropylene components.

Typical applications: Dry heat sterilization, autoclaving of liquid media, storage of corrosive reagents and sampling.

Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	h (mm)	Pack Unit
<b>With high temperature screw cap and ETFE pouring ring</b>					
21 801 63 16	2 000	45	136	265	10
21 801 73 12	5 000	45	182	335	1
21 801 86 17	10 000	45	227	415	1



> Find your nearest **distributor** on our global network:  
[www.DWK-LifeSciences.com/DURAN/distributors](http://www.DWK-LifeSciences.com/DURAN/distributors)



DURAN® PROTECT LABORATORY BOTTLE



with DIN thread, plastic coated



With easy-to-read scale. In fired-on, highly durable white ceramic. Service temperature limit of the PU plastic coating: -30 °C to +135 °C. The coating provides scratch, leak\* and splinter\* protection and is ideally suited to both the transport and storage of hazardous media or valuable samples. UV protection up to approx. 380nm wavelength. High transparency. Suitable for microwaving. (\* only applies to bottles 5000 mL and less)

Typical applications: storage, transport and safe handling of hazardous or valuable substances.

Cat. No.	Capacity (mL)	DIN Thread (GL)	d (OD) (mm)	h (mm)	Pack Unit
<b>with screw cap and pouring ring from PP (blue)</b>					
21 805 63 52	2 000	45	136	260	10
21 805 73 57	5 000	45	182	330	1
<b>without screw cap and pouring ring</b>					
21 805 63 03	2 000	45	136	260	10
21 805 73 08	5 000	45	182	330	1
21 805 86 04	10 000	45	228	410	1
21 805 91 03	20 000	45	289	505	1



> Find your nearest distributor on our global network:  
[www.DWK-LifeSciences.com/DURAN/distributors](http://www.DWK-LifeSciences.com/DURAN/distributors)

